

Righton's Lightweight Test Frames

General specifications (All models can hold three lenses at the front and one at the back).

For men



White PD: 54 mm–70 mm (DXE30007)

Gray PD: 54 mm–70 mm (DXE30006)


For women and children



Black & White PD: 48 mm–60 mm (DXE30008)


Gray PD: 48 mm–60 mm (DXE30005)

Features

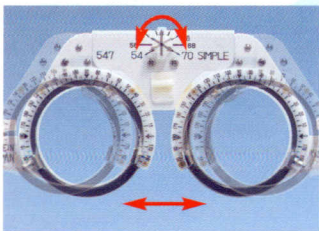


Adjustable pupil distance

Turn the dial to select pupil distances from 54 mm to 70 mm or 48 mm to 60 mm, depending on model.

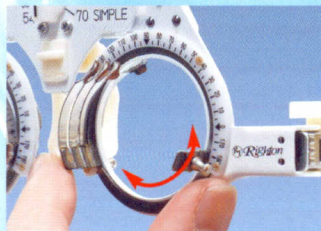


Axis of cylinder

Rotate the frame ring by five-axis-step clicks to adjust the axis.

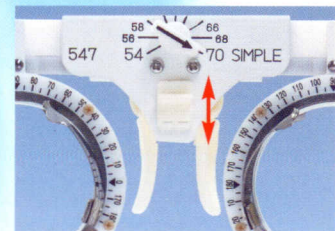


Progressive lens support

Move the switch down to provide more accurate testing when using progressive lenses.

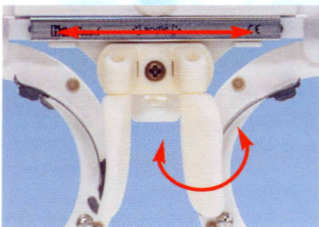


Adjustable nose pad

For comfort, move the nose pad end piece to either one of two positions.

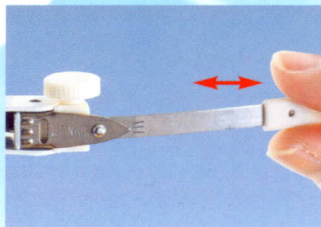

Extending temple tips

To adjust temple tips, hold each hinge and temple tip and pull in opposite directions. Each temple tip can be moved a maximum of 10 mm, sufficient for most head sizes (30 mm on the 5470 model and 20 mm on the 4860 model).

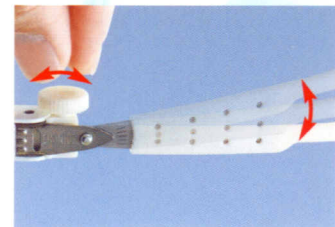

Temple inclination

Simply rotate the small dial on the temple to adjust vertically by a maximum of 14 degrees.